

# **Elmasonic Easy**

# **Product Information and Prices**





### **Contents**

Elmasonic Easy	3
Order numbers and list prices	
Product features and functions	
Technical Data - Mechanics	. 5
Technical Data - Electronics and noise levels	. 6
Dimension definition	.7
Elmasonic Easy R	9
Order numbers and list prices	
Product features and functions	
Technical Data - Mechanics	
Technical Data - Electronics and noise levels	11
Dimension definition	12



### **Elmasonic Easy**

Ultrasonic cleaning devices



#### Order numbers and list prices

Elmasonic Easy	10*	10H	20*	20H	30*	30H	40H	60H	100H	120H	180H	300H
Order number 220 - 240 V~	1118881	1118883	1118885	1118887	1118889	1117697	1118892	1118894	1118896	1118898	1118900	1118902
Order number 115 - 120 V~	1118880	4440000	1118884	4440000	1118888	1118890	1118891	4440000	1118895	1118897	1118899	1118901
Order number 100 V~		1118882		1118886		1126018		1118893				
List price	305,00€	370,00€	365,00 €	440,00€	420,00€	470,00€	695,00€	810,00€	970,00€	1.260,00 €	1.370,00 €	1.720,00€

#### Scope of delivery

- Ultrasonic unit
- Cover or, turned over, usable as a drip tray for a basket
- Hose connector for drain (from Elmasonic Easy 60)
- Power cable with rubber connector (2 m length)
- Manual

\*Unit without heating

All prices ex works, excluding packaging.



### **Elmasonic Easy**

#### **Product features and functions**

- 37 kHz ultrasonic frequency (control range -2,5 / +5,5 kHz)
- 2 ultrasonic modes for need-based sound field distribution in the ultrasonic bath
- dynamic mode, automatically set: Serves more intensive cleaning with higher cleaning power. Sweep and pulse alternately for homogeneous and powerful cleaning.
- eco mode, can be switched on: Serves gentle cleaning with quieter operation. Saves energy and increases the life of the tub.
- Electronic selector switches on wipeable membrane keypad and LED display of all operating states
- Pre-selection of cleaning time: cleaning time adjustable between 1 minute and 30 minutes or continuous operation adjustable (ultrasonic continuous operation max. 6 hours)
- Acoustic and optical signal at the end of cleaning
- Electronic temperature control adjustable in steps of 5 °C from 30 80 °C for devices with heating (temperature tolerance -5 / +8 °C)
- Safety switch-off with wake-up function
- after 8 hours of inactivity
- in case of a brief interruption of the power supply
- Automatic shutdown of ultrasonic operation at a bath temperature > 90 °C
- Splash-proof housing and cleaning tray made of stainless steel
- Side stopcock for emptying the cleaning tank via the drain at the rear of the unit (from Elmasonic Easy 60H)
- 2-year warranty from the initial operation, but no longer than 3 years from the production date



### **Elmasonic Easy**

#### **Technical Data - Mechanics**

Elmasonic Easy	10	10H	20	20H	30	30H	40H	60H
Total capacity V (I)	0,9		1,6		2,7		3,9	5,9
Recom. operating capacity V <sub>a</sub> (I)	0	,7	1,2		1,6		2,9	4,3
Unit ext. dimensions W / D / H (mm)	235 / 13	35 / 200	215 / 1	215 / 190 / 230 305		190 / 230 305 / 190 / 280		370 / 225 / 280
Opening dimensions W <sub>2</sub> / D <sub>2</sub> (mm)	190	190 / 85		150 / 135 235 / 135		235 / 135	300 / 150	
Tank int. dimensions W <sub>1</sub> / D <sub>1</sub> / H <sub>1</sub> (mm)	185 / 8	80 / 55	140 / 125 / 95		220 / 120 / 95		215 / 115 / 145	280 / 130 / 145
Weight incl. cover (kg)	2	,0	2,1		3,3		3,3 4,0	
Drain hose - ID (mm)						12,0		
Basket int. dimensions W <sub>3</sub> / D <sub>3</sub> / H <sub>3</sub> (mm)	155 / 7	70 / 35	110 / 1	105 / 50	195 / 1	05 / 50	195 / 105 / 95	255 / 105 / 95
Total basket loading (kg)	1,	,0	1	,0	1,	0	3,0	5,0

Elmasonic Easy	100H	120H	180H	300H
Total capacity V (I)	9,5	12,9	17,8	27,5
Recom. operating capacity V <sub>a</sub> (I)	6,8	10,0	14,2	22,0
Unit ext. dimensions W / D / H (mm)	380 / 315 / 280	380 / 315 / 330	405 / 375 / 330	575 / 375 / 330
Opening dimensions W <sub>2</sub> / D <sub>2</sub> (mm)	295 / 235	295 / 235	325 / 295	500 / 295
Tank int. dimensions W <sub>1</sub> / D <sub>1</sub> / H <sub>1</sub> (mm)	280 / 220 / 145	275 / 215 / 195	300 / 270 / 195	475 / 270 / 195
Weight incl. cover (kg)	5,9	7,5	8,5	11,0
Drain hose - ID (mm)		12	2,0	
Basket int. dimensions W <sub>3</sub> / D <sub>3</sub> / H <sub>3</sub> (mm)	250 / 195 / 90	250 / 195 / 135	280 / 255 / 135	450 / 255 / 135
Total basket loading (kg)	6,0	7,0	8,0	10,0



### **Elmasonic Easy**

#### **Technical Data - Electronics and noise levels**

Elmasonic Easy	10	10H	20	20H	30	30H	40H	60H	100H	120H	180H	300H
Mains frequency (Hz)		50 / 60										
Power consumption total (W)	30	90	35	155	80	280	320	55	50	10	000	1500
Ultrasonic power effective ± 10 % (W)	3	30	3	35	8	0	120	150 20		00	300	
Ultrasonic peak power max.** (W)	24	40	28	80	32	20	480	60	00	8	00	1200
Heating power (W)	0	60	0	120	0	200	200	40	00	8	00	1200
Temperature (°C)	25 - 80											
Protection class	IP 20											





### **Elmasonic Easy**

#### **Dimension definition**

Unit outer dimensions


Unit sizes 10 - 20

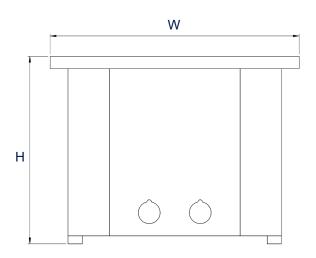


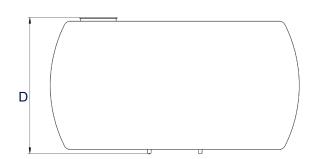

#### Unit outer dimensions


Unit sizes 30 - 300

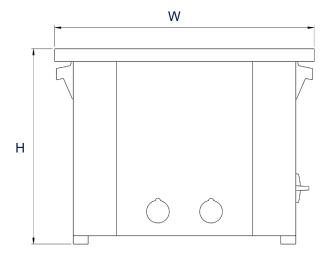


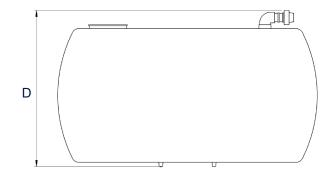

#### Tank internal dimensions

 $V_{\rm a}$ : recommended operating capacity of the cleaning tank (The actual operating capacity depends on the load.)

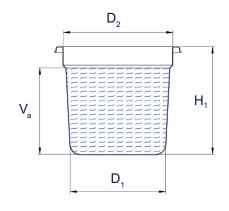

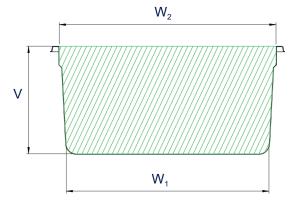


The drawings are meant as a visualization aid and are not true to scale.




### **Elmasonic Easy**

#### **Dimension definition**

Basket internal dimensions

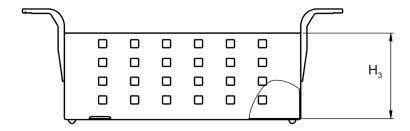

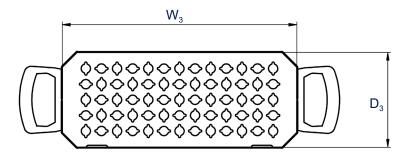



### **Elmasonic Easy R**

Ultrasonic cleaning devices



#### Order numbers and list prices

Elmasonic Easy	50 R	350 R
Order number 220 - 240 V~	1120119	1120120
Order number 100 - 120 V~	1120118	
List price	955,00 €	3.190,00 €

#### Scope of delivery

- Ultrasonic unit
- Cover or, turned over, usable as a drip tray for a basket (50 R)
- Stainless steel lid (350 R)
- Hose connector for drain (350 R)
- Power cable with rubber connector (2 m length)
- Manual

All prices ex works, excluding packaging.



### Elmasonic Easy R

#### **Product features and functions**

- 37 kHz ultrasonic frequency (control range -2,5 / +5,5 kHz)
- 4 ultrasonic modes for demand-oriented sound field distribution and sample preparation in the ultrasonic bath
- general cleaning sweep mode: General cleaning tasks. Ensures a more homogeneous distribution of cavitation (mechanical cleaning action) in the bath.
- sieve cleaning mode: Increased ultrasonic intensity for the most stubborn soiling, saves time. For optimal cleaning of test sieves.
- sample preparation mode: Sweep and pulse alternately for homogeneous and strong cleaning. For sample preparation such as mixing, dissolving, dispersing, cleaning, etc.
- degas mode: more effective, faster and more reliable. Fast degassing of HPLC solvents
- Electronic selector switches on wipeable membrane keypad and LED display of all operating states
- Pre-selection of cleaning time: cleaning time adjustable between 1 minute and 30 minutes or continuous operation adjustable (ultrasonic continuous operation max. 6 hours)
- Acoustic and optical signal at the end of cleaning
- Safety switch-off with wake-up function
- after 8 hours of inactivity
- in case of a brief interruption of the power supply
- Automatic shutdown of ultrasonic operation at a bath temperature > 90 °C
- Splash-proof housing and cleaning tray made of stainless steel
- Stopcock for emptying the cleaning tank via the drain at the rear of the unit (350 R)
- Different fill level possible for different container sizes. Saving of water and medium (50 R)
- 2-year warranty from the initial operation, but no longer than 3 years from the production date

09 15 Subject to technical and visual modificat



### **Elmasonic Easy R**

#### **Technical Data - Mechanics**

Elmasonic Easy	50 R	350 R
Total capacity V (I)	6,0	69,0
Recom. operating capacity V <sub>a</sub> (I)	3,0 - 5,0	59,0
Unit ext. dimensions W / D / H (mm)	345 / 295 / 345	675 / 795 / 505
Tank int. dimensions Dia / H <sub>1</sub> (mm)	Ø 240 / 140	Ø 595 / 245
Weight incl. cover (kg)	5,0	21,0
Drain hose - ID (mm)		12,0
Basket int. dimensions W <sub>2</sub> / H <sub>2</sub> (mm)	Ø 220	Ø 545
Total basket loading (kg)	4,0	10,0

#### **Technical Data - Electronics and noise levels**

Elmasonic Easy	50 R	350 R				
Mains voltage (V~)	115 - 120 / 220 - 240	220 - 240				
Mains frequency (Hz)	50 / 60	50				
Power consumption total (W)	150	650				
Ultrasonic power effective ±10% (W)	150	650				
Ultrasonic peak power max.* (W)	600	2600				
Protection class	IP 20					
Sound pressure level (LpAU)**	< 80 dB					

<sup>\*</sup> The signal form of the wave results in a factor 4 for the ultrasonic peak max., depending on the modulation of the wave.

<sup>\*\*</sup> Measured sound pressure level with basket and cover at 1 m distance.





### **Elmasonic Easy R**

#### **Dimension definition**


Unit outer dimensions


Unit size 50R

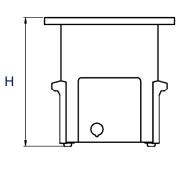

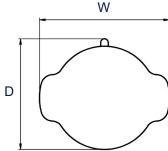

Unit size 350R

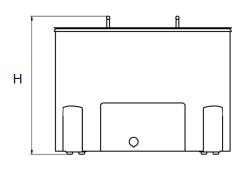

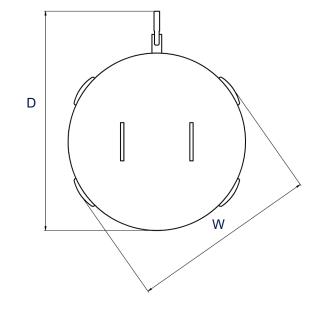



#### Tank internal dimensions

V<sub>a</sub>: recommended operating capacity of the cleaning tank (The actual operating capacity depends on the load.)



#### Basket internal dimensions

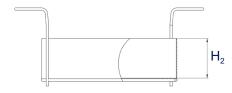

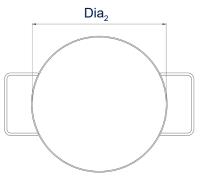

The drawings are meant as a visualization aid and are not true to scale.